



3003 Imaging Table



*Multi-Directional
Top Slide for
Optimal Full Body
Image Coverage*

An Imaging Table Built to Skytron's Quality Engineering & Design Standards

Skytron's Impulse 3003 Imaging Table is built with the same high quality, precision design and durability standards as our Skytron surgical tables. Impulse 3003 is equipped with an 83 in. image-friendly, extended carbon fiber top offering 2 Speed, Multi-directional Top Slide for optimal Full Body Imaging Coverage.




Joy Stick Controller

- 2 Speed, Multi-directional Carbon Fiber Top
- 500 lb. Weight Capacity
- Narrow Table Base Profile for Ease of Patient Access
- Automatic 4 point, Self Leveling Brakes for Optimal Stability
- Durable Back-lit, One Touch Pendant Control & Multi-directional Joy Stick Controller
- Each Controller is Operable from Any Position Around the Table's Periphery

Impulse 3003 Imaging Table with 2 Speed, Multi-directional Top Slide

Impulse 3003 provides expanded patient positioning options with power assisted 2 Speed, Multi-directional Top Slide, eliminating any need to move the C-arm.

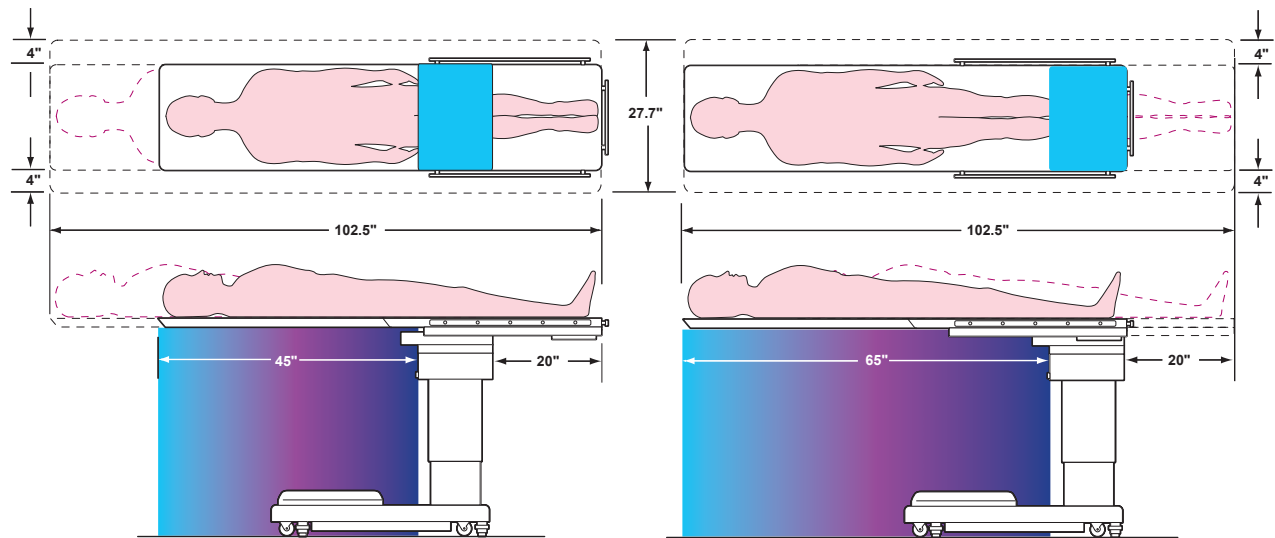

Imaging Area:

20 in. Longitudinal and 8 in. Lateral movement, providing an expanded and mobile imaging window of 27.7 in. X 65 in., along with physician controlled, 2 speed, multi-directional top slide.

- 83 in. Carbon Fiber Imaging Surface
- 8 in. Lateral Slide, 20 in. Longitudinal Slide
- Unobstructed Fluoroscopic Imaging Area of 65 in.
- Multi-directional Joy Stick Controller with Reverse Mode, High/Low Speed and Auto Centering Buttons



42 in. Max
Table Height



28 in. Min
Table Height



20°
Trendelenburg



20° Reverse
Trendelenburg



20 in. Longitudinal
Top Slide



8 in. Lateral
Top Slide



20° Lateral Tilt
Right & Left

3003 Imaging Table Specifications

Technical Specifications:

Table Top

Length	83 in.	2100mm
Width	19.7 in.	500mm
Width with side rails	22 in.	560mm

Adjustments via hand control:

Height	28-42 in.	710-1065mm
Lateral tilt (left/right)	20°	
Trendelenburg	20°	
Reverse Trendelenburg	20°	
Return to level (return to horizontal position from articulated positions)		
Base braking (lock/unlock)		

Joy Stick Control:

Top button on Joy Stick Control must be depressed and held while moving Joy Stick to operate

Orientation Button:	Reverse	Normal
Center Button:	Moves table top to center position	
Speed:	High	Low

Weight of operating table 500 lb. (227kg)

Max. permissible load 500 lb. (227kg) Lift
500 lb. (227kg) Articulation



Full function, electrohydraulic surgical table including 500 lb. (227kg) lift and 500 lb. (227kg) articulation capacity. Table is equipped with a four point, self leveling brake system to assure stability and safety. A manual emergency brake release system is provided to allow moving the table in the event of a power loss. Single touch pendant and multi-directional joystick controllers, connected to the table base by coiled cords, permits 2 speed, multi-directional top slide and functional operation from any point around the table's perimeter.



5085 Corporate Exchange Blvd. SE • Grand Rapids, MI 49512
1-800-SKYTRON (759-8766) • email: sales@skytron.us

Local Representative

3/2008