




# 3502 EZ SLIDE

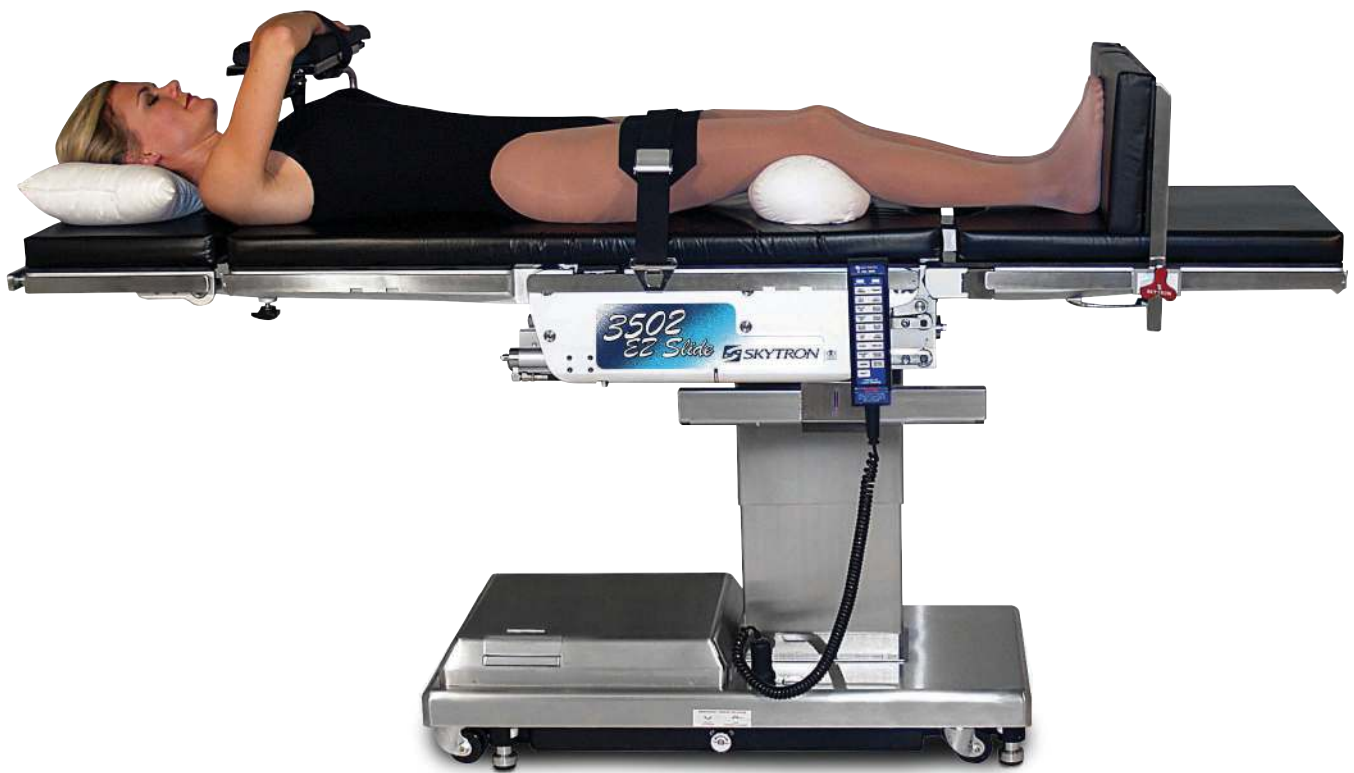

 **SKYTRON**<sup>®</sup>  
TRUSTED PARTNER | PROVEN SOLUTIONS

# Simplify full body imaging and patient positioning.

Featuring convenient top slide and the tallest table height in our portfolio, the 3502 EZ Slide is designed with the surgical team in mind. Imaging and surgical access is easy with 21" of top slide, allowing for movement of the patient through the imaging window with just the click of a button. Also featuring our tallest table height, the 3502 EZ Slide provides comfortable waist-height access, minimizing the physical strain placed on the surgical team.

## **IMAGING EFFICIENCY**

With 21" of top slide, the 3502 EZ Slide makes full body imaging convenient. A simple click of a button effortlessly moves the patient toward the C-arm, eliminating the need for risky repositioning. The table also features a slim profile base that stays clear of mobile imaging equipment.

## **SURGICAL POSITIONING FLEXIBILITY**

Skytron makes positioning and setup a breeze with a removable leg section that can quickly integrate with accessories such as the optional 40" carbon fiber extension. Sturdy, standardized side rails accept most table accessories and preprogrammed positioning options allow for one-touch maneuvers at any point throughout the procedure.

## **IMPROVED PATIENT SAFETY**

The 3502 EZ Slide takes the confusion out of weight limits. Skytron can provide top slide at full weight capacity, which means the 3502 EZ Slide will support articulation capacities up to 600 pounds – even while sliding the patient through the imaging window.

**COMFORTABLE  
SURGICAL  
ACCESS**

**FULL  
ARTICULATION  
CAPACITY**

**ONE-TOUCH  
IMAGING  
ACCESS**



## 3502 EZ Slide

### 21" TOP SLIDE, GENERAL PURPOSE TABLE

Providing 21" of top slide, the 3502 EZ Slide takes the work out of imaging setup. At the click of a button, the patient is moved through the expansive imaging window.



### DIFFERENTIATING FEATURES

- 21" top slide
- 7<sup>5</sup>/<sub>8</sub>" slide to head and 13<sup>1</sup>/<sub>8</sub>" slide to foot
- 700 lb lift and 600 lb articulation weight capacity
- Removable leg section for quick accessory attachment and positioning
- Durable stainless steel base is easy to clean
- 26"- 44<sup>7</sup>/<sub>8</sub>" adjustable table height
- 20° lateral tilt
- 25° Trendelenburg and 20° reverse Trendelenburg





## Convenient and Comfortable Imaging Access

Skytron makes top slide easy for patients of diverse sizes. The 3502 EZ Slide will offer 21" of top slide at its full articulation weight capacity. Also designed to give comfortable, waist-height access to taller surgeons, the table offers our maximum table height of 44 <sup>7</sup>/<sub>8</sub>". As the table top slides toward the foot end, the head section can be conveniently detached and applied as a foot extension. This thoughtful design is complemented with a radiolucent top and a removable leg section that allows for quick attachment of the 40" carbon fiber extension so that any operating room can achieve improved imaging without the worry of reduced weight limits.

The 3502 EZ Slide also brings confidence to the operating room. Recognized for their reliable performance, Skytron surgical tables are consistently rated as having a lower cost of ownership.



26" Minimum  
Table Height



44 <sup>7</sup>/<sub>8</sub>" Maximum  
Table Height



21" Total Slide



40° Flex



90° Reflex



98° Leg Down



One-touch positioning functions include Trendelenburg, tilt and more



Battery life/diagnostic indicator



Easy to use clip style power cord



Back-up controls for emergency situations



Optional 40" carbon fiber leg section for upper and lower body imaging



Self-leveling brakes for stability on uneven floors

## ONE-TOUCH PENDANT CONTROL

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Audible alarm for kidney bridge, back up articulation and top slide obstruction
- Anatomically correct return-to-level
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



90° Back Up



40° Back Down



25° Trendelenburg




20° Reverse Trendelenburg



20° Tilt - Left



20° Tilt - Right

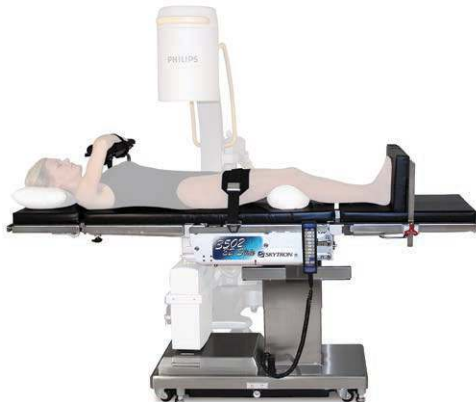

UPPER BODY IMAGING



LOWER BODY IMAGING



LAP NISSEN



SPINAL

## Positioning Guide



OPHTHALMIC/ENT



NEPHRECTOMY



NEURO (NECK)



NEURO (LUMBAR)



SHOULDER ARTHROSCOPY



UROLOGY



LOWER EXTREMITY



BARIATRIC



## 3502 EZ SLIDE

### SPECIFICATIONS

Length w/o head rest	69 3/8"	1,750mm
Length with head rest	82"	2,083mm
Width (back/seat plate)	19 5/8"	500mm
Width with side rails	21 7/8"	555mm

#### Adjustments via pendant control:

Height	26" - 44 7/8"	660-1,140mm
Lateral tilt (left/right)	20°	
Trendelenburg	25°	
Reverse Trendelenburg	20°	
Power Kidney Bridge	5 7/8"	150mm Lift
Back up/down	90°/40°	
Leg up/down	0°/98°	
Return to level (anatomically correct return to horizontal position from articulated positions)		
Top slide/head	7 5/8"	194mm
Top slide/foot	13 1/8"	335mm
Flex/Reflex	40°/90°	
Base braking (lock/unlock)		

#### Adjustments via manual control:

Head Section	60° up/90° down	
Max. permissible load	700 lb (315kg)	Lift
	600 lb (270kg)	Articulation

Full function, electrohydraulic battery powered surgical table features 21" longitudinal top slide, 700 lb lift and 600 lb articulation capacity. Table is equipped with a removable 15' power cord for battery charging and removable leg section for accessory attachment. Safety features include four conductive casters, four-point braking system with push button operation and manual override. Handheld, backlit pendant control includes steel housing and provides extreme durability. Pendant permits control at any point around the table's perimeter.

Printed with agra-based inks, water based coatings on FSC®-certified paper.



5085 Corporate Exchange Blvd. SE  
Grand Rapids, MI 49512  
1.800.SKYTRON (759-8766)  
email: [info@skytron.us](mailto:info@skytron.us)  
[www.skytron.us](http://www.skytron.us)

**About Skytron** Skytron is a privately held company specializing in capital equipment that drives efficiency in today's healthcare facilities. Our surgical tables, surgical & exam lighting, booms, video integration systems and other products help hospitals around the world work better. To learn more, visit our website at [www.skytron.us](http://www.skytron.us).