



# 3602 ULTRASLIDE



 **SKYTRON**<sup>®</sup>  
TRUSTED PARTNER | PROVEN SOLUTIONS

# Delivering 23" of top slide without any compromises.

Pairing the industry's longest top slide with substantial weight capacities, imaging access is only part of the 3602 UltraSlide advantage. Combine the table's top slide capabilities with an 800 pound articulation capacity, and the result is the best of both worlds: imaging and power. The table's impressive 23" of top slide accommodates the industry's most generous articulation capacity and provides superior imaging for patients of diverse sizes.

## **IMAGING EFFICIENCY**

The 3602 UltraSlide makes full-body imaging easier. At the touch of a button, patients can be effortlessly moved through the imaging window. The low-profile base stays clear of the imaging equipment, providing hassle-free C-arm access.

## **SURGICAL POSITIONING FLEXIBILITY**

The removable back and leg section allow for easy and efficient accessory attachment. Enhance imaging capabilities with the optional 40" carbon fiber extension or attach the beach chair positioner for shoulder surgeries. With the back section removed, the beach chair can be fully integrated and powered by the table, providing one-touch setup.

## **IMPROVED PATIENT SAFETY**

Capable of 23" top slide at a full articulation capacity of 800 pounds, the 3602 UltraSlide will remain secure for patients of larger sizes. Its extensive capabilities not only make patient imaging simple, but also eliminate the need for risky repositioning

**MARKET-  
LEADING  
TOP SLIDE**

**POWERFUL  
ARTICULATION  
CAPACITY**


**ONE-TOUCH  
IMAGING  
ACCESS**



## 3602 UltraSlide

### 23" TOP SLIDE, GENERAL PURPOSE TABLE

Offering the industry's maximum top slide, the 3602 UltraSlide takes the work out of imaging. At the click of a button, the patient can be moved throughout the expansive imaging window.

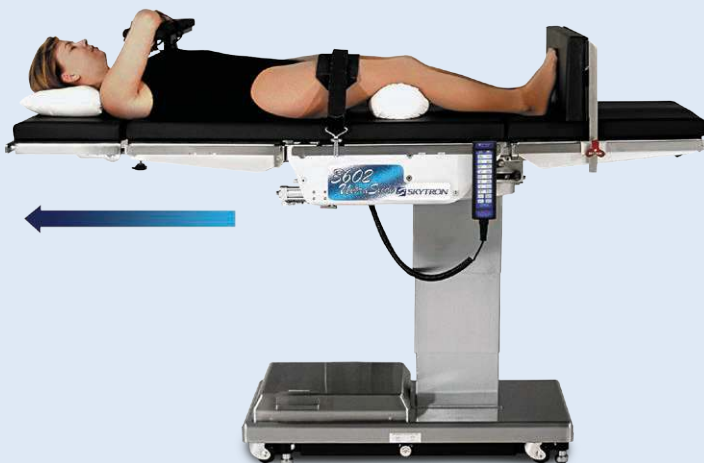

### DIFFERENTIATING FEATURES

- 23" top slide
- 9 1/2" slide to head and 13 1/4" slide to foot
- 1,000 lb lift and 800 lb articulation weight capacity
- Removable back and leg sections for quick accessory attachment and positioning
- Power beach chair positioner for one-touch set-up
- Durable stainless steel base is easy to clean
- 24 3/8"-43 1/4" adjustable table height
- 30° lateral tilt
- 30° Trendelenburg and reverse Trendelenburg





## More Top Slide and Power in One Table

Able to support a 1,000 pound lift and 800 pound articulation weight capacity, the 3602 UltraSlide provides 23" of top slide with greater power. Its superior design supports the full articulation capacity even during 23" of top slide – something the competition cannot do. As the table top slides toward the foot end, the head section can be conveniently detached and applied as a foot extension. With removable back and leg sections, the table offers quick attachment of accessories like the 40" carbon fiber extension and power beach chair positioner. Pre-programmed positioning options make setup even easier, allowing for one-touch maneuvers at any point throughout the procedure.

The 3602 UltraSlide also brings confidence to the operating room. Recognized for their reliable performance, Skytron surgical tables are consistently rated as having a lower cost of ownership.



24 3/8" Minimum  
Table Height



43 1/4" Maximum  
Table Height



40° Flex



90° Reflex



15° Leg Up



100° Leg Down



One-touch positioning functions include return-to-level and beach chair



Battery life/diagnostic indicator



Easy to use clip style power cord



Back-up controls for emergency situations



Optional 40" carbon fiber leg section for upper and lower body imaging



Self-leveling brakes for stability on uneven floors

## ONE-TOUCH PENDANT CONTROL

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Audible alarm for kidney bridge, back up articulation and top slide obstruction
- Anatomically correct return-to-level
- Anatomically correct beach chair position
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



90° Back Up



40° Back Down



30° Trendelenburg



30° Reverse Trendelenburg



30° Tilt - Left



30° Tilt - Right

Positioning Guide



UPPER BODY IMAGING



LOWER BODY IMAGING



UPPER BODY IMAGING WITH  
40" CARBON FIBER EXTENSION



LOWER BODY IMAGING WITH  
40" CARBON FIBER EXTENSION



OPHTHALMIC/ENT



KIDNEY/THORACIC

Bariatric Guide



ABDOMINAL



OPHTHALMIC





NEURO (NECK)



NEURO (LUMBAR)



SHOULDER ARTHROSCOPY



LAP NISSEN



UROLOGY



LUMBAR



GASTRO/INTESTINAL



CYSTO/GYN/ENDUROLOGY



## 3602 ULTRASLIDE

### SPECIFICATIONS

Length w/o head rest	69 1/2"	1,765mm
Length with head rest	82 1/8"	2,086mm
Width	19 5/8"	500mm
Width with side rails	21 7/8"	555mm

#### Adjustments via pendant control:

Height	24 3/8" to 43 1/4"	620-1,100mm
Lateral tilt (left/right)	30°	
Trendelenburg	30°	
Reverse Trendelenburg	30°	
Power Kidney Bridge	5 3/4"	146mm Lift
Back up/down	90°/40°	
Leg up/down	15/100°	
Return to level (anatomically correct return to horizontal position from articulated positions)		
Top slide/head	9 1/2"	245mm
Top slide/foot	13 1/4"	335mm
Flex/Reflex	40°/90°	
Base braking (lock/unlock)		
Beach Chair		

#### Adjustments via manual control:

Head Section	60° up/90° down	
Max. permissible load	1,000 lb (450kg)	Lift
	800 lb (360kg)	Articulation

Full function, electrohydraulic battery powered surgical table features 23" longitudinal top slide, 1,000 lb lift and 800 lb articulation capacity. Table is equipped with a removable 15' power cord for battery charging and removable back and leg sections for accessory attachment. Safety features include four conductive casters, four-point braking system with push button operation and manual override. Handheld, backlit pendant control includes steel housing and provides extreme durability. Pendant permits control at any point around the table's perimeter.

Printed with agra-based inks, water based coatings on FSC®-certified paper.



5085 Corporate Exchange Blvd. SE  
 Grand Rapids, MI 49512  
 1.800.SKYTRON (759-8766)  
 email: [info@skytron.us](mailto:info@skytron.us)  
[www.skytron.us](http://www.skytron.us)

**About Skytron** Skytron is a privately held company specializing in capital equipment that drives efficiency in today's healthcare facilities. Our surgical tables, surgical & exam lighting, booms, video integration systems and other products help hospitals around the world work better. To learn more, visit our website at [www.skytron.us](http://www.skytron.us).