



6702 HERCULES



 **SKYTRON**[®]
TRUSTED PARTNER | PROVEN SOLUTIONS

Giving operating rooms the ability to do more.

Combining impressive 1,200 pound lift capacity with exclusive table top rotation, the 6702 Hercules exceeds all standards in general purpose surgical tables. Skytron offers 210° table top rotation – making imaging and patient positioning easy and safe. From pediatric to bariatric, the 6702 Hercules gives any operating room the versatility to handle more procedures and specialties.

IMAGING EFFICIENCY

At the turn of a handle, the 6702 Hercules allows superior imaging access at both the head and foot end of the table. The table can quickly rotate toward the C-arm and the low-profile base stays clear of the imaging equipment. Skytron takes the work out of imaging setup.

SURGICAL POSITIONING FLEXIBILITY

A remarkable 1,200 pound lift and 1,000 pound articulation capacity makes the 6702 Hercules ready for patients of diverse sizes. It also features the lowest table height in our portfolio. The removable back and leg sections allow seamless integration with accessories such as the power beach chair positioner – providing automatic, one-touch setup for shoulder surgery.

IMPROVED PATIENT SAFETY

With a table top that rotates independent of the base, the 6702 Hercules allows the brakes to remain firmly in place while moving the patient and maximizing surgical access. This unique feature sets the 6702 Hercules apart and makes any operating room safe.

**EFFICIENT
TOP
ROTATION**

**GENEROUS
WEIGHT
CAPACITY**

**INCREASED
SURGICAL
SPECIALTIES**



6702 Hercules

210° TOP ROTATION, GENERAL PURPOSE TABLE

Skytron offers true top rotation, allowing the surgeon or C-arm to remain at the foot end of the table where the base does not interfere. And with the removable back and leg section, the power beach chair can be attached for the most comfortable shoulder surgery setup available.



DIFFERENTIATING FEATURES

- 210° top rotation
- 1,200 lb lift and 1,000 lb articulation weight capacity
- Removable back and leg sections for quick accessory attachment and positioning
- Power beach chair positioner for one-touch set-up
- Durable stainless steel base is easy to clean
- 23 5/8"-41 1/2" adjustable table height
- 30° lateral tilt
- 30° Trendelenburg and reverse Trendelenburg





One Table, More Specialties

The 6702 Hercules allows 210° top rotation at the turn of a handle, making procedure and imaging setup easier than ever. Seamlessly rotate the table toward the C-arm or transition the patient between the surgeon and anesthesiologist during head-end procedures. Featuring 1,200 pound lift and 1,000 pound articulation capacity combined with our lowest table height, the 6702 Hercules can also support bariatric patients with comfortable surgical access. Its radiolucent top permits sharp image amplification while the sturdy side rails will accept a full range of complementary accessories. From its ease of setup to its impressive power, the 6702 Hercules helps any operating room manage a greater variety of procedures with only one surgical table!

The 6702 Hercules also brings confidence to the operating room. Recognized for their reliable performance, Skytron surgical tables are consistently rated as having a lower cost of ownership.



23 5/8" Minimum
Table Height



41 1/2" Maximum
Table Height



40° Flex



90° Reflex



20° Leg Up



100° Leg Down



One-touch positioning functions include return-to-level and beach chair



Battery life/diagnostic indicator



Easy to use clip style power cord



Optional side extensions widen the table top to 30"



Optional 40" carbon fiber leg section for upper and lower body imaging



Self-leveling brakes for stability on uneven floors

ONE-TOUCH PENDANT CONTROL

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Anatomically correct return-to-level
- Anatomically correct beach chair position
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



90° Back Up



40° Back Down



30° Trendelenburg



30° Reverse Trendelenburg



30° Tilt - Left



30° Tilt - Right



UPPER BODY IMAGING



LOWER BODY IMAGING



LAP NISSEN



UROLOGY



OPHTHALMIC/ENT



NEPHRECTOMY



ABDOMINAL




NEURO (NECK)



NEURO (LUMBAR)



SHOULDER ARTHROSCOPY

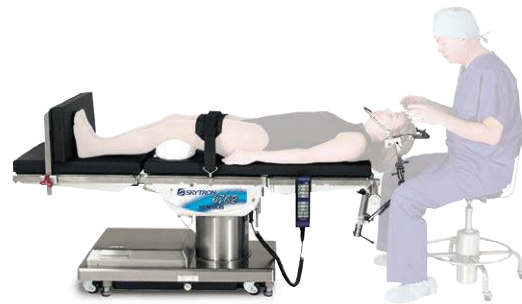

NEURO (HEAD)



HIP REVISION



CARDIOVASCULAR



GALL BLADDER



CYSTO/GYN



6702 HERCULES

SPECIFICATIONS

Length w/o head rest	66 3/8"	1,665mm
Length with head rest	79"	2,006mm
Width (back / seat plate)	19 5/8"	500mm
Width with side rails	21 7/8"	555mm

Adjustments via pendant control:

Height	23 5/8 to 41 1/2"	600-1,053mm
Lateral tilt (left/right)	30°	
Trendelenburg	30°	
Reverse Trendelenburg	30°	
Power Kidney Bridge	5 7/8" Lift	150mm Lift
Back up/down	90°/40°	
Leg up/down	20°/100°	
Return to level (anatomically correct return to horizontal position from articulated positions)		
Beach Chair		
Flex/Reflex	40°/90°	
Base braking (lock/unlock)		

Adjustments via manual control:

Top Rotation	210°
Head Section	60° up/90° down

Max. permissible load	1,200 lb (545kg)	Lift
	1,000 lb (455kg)	Articulation

Full function, electrohydraulic battery powered surgical table features 210° top rotation, 1,200 lb lift and 1,000 lb articulation capacity. Table is equipped with a removable 15' power cord for battery charging, removable back and leg section for accessory attachment, and it exhibits the lowest table height in our portfolio. Safety features include four conductive casters, four-point braking system with push button operation and manual override. Handheld, backlit pendant control includes steel housing and provides control at any point around the table's perimeter.

Printed with agra-based inks, water based coatings on FSC®-certified paper.



5085 Corporate Exchange Blvd. SE
 Grand Rapids, MI 49512
 1.800.SKYTRON (759-8766)
 email: info@skytron.us
 www.skytron.us

About Skytron Skytron is a privately held company specializing in capital equipment that drives efficiency in today's healthcare facilities. Our surgical tables, surgical & exam lighting, booms, video integration systems and other products help hospitals around the world work better. To learn more, visit our website at www.skytron.us.