



Stainless Warming Cabinets



J2 Series

Skytron Warming Cabinets

Skytron's J2 Series Warming Cabinets are designed to maintain clinically ideal temperatures for solution bottles, solution bags and blankets. Six models are available to provide optimum clinical storage requirements.

Free Standing



Split Glass Door
SS2207-J2G

Counter Top	
SS2201-J2 (solid door)	Twenty-eight (28) 1500 ml solution bottles or eleven (11) blankets
SS2201-J2G (glass door)	

Free Standing		
SS2205-J2	Single compartment (single solid door)	Twenty-eight (28) 1500 ml solution bottles or eleven (11) blankets per shelf, including separate upper compartment.
SS2205-J2G	Single compartment (single glass door)	
SS2207-J2	Two compartment (split solid doors)	
SS2207-J2G	Two compartment (split glass doors)	

Model	Watts	Amps	Weight
SS2201-J2	240	2.0	179 lbs.
SS2201-J2G	240	2.0	206 lbs.
SS2205-J2	720	6.5	340 lbs.
SS2205-J2G	720	6.5	391 lbs.
SS2207-J2	1056	8.8	404 lbs.
SS2207-J2G	1056	8.8	464 lbs.

Counter Tops



Glass Door
SS2201-J2G



Solid Door
SS2201-J2

Standards:

- ETL Listed for US and Canada
- Seismic Anchorage Calculations with OSHPOD Approved R Numbers

Free Standing



Split Solid Doors
SS2207-J2

Single Glass Door
SS2205-J2G

Single Solid Door
SS2205-J2

Construction:

- Double wall stainless steel construction with number 4 finish
- Field reversible left or right hinged doors
- Optional double paned, tempered glass doors
- Heating compartments fully insulated with circulating fan
- Two or three adjustable stainless shelves (depending upon model)

Operating & Safety Controls:

- Digital temperature controls for individual compartments
- Recirculating fan in each compartment
- Adjustable temperature range
- Set point temperature +/- 5° F
- Independent power switch / compartment
- Audible and visual over-temperature indicator
- Automatic shut off (temperature exceeding 180° F)

Innovating Efficient Solutions for Healthcare Since 1972



Hybrid Rooms




Surgical, Procedure
and Exam Lights



Tables



Financing Options

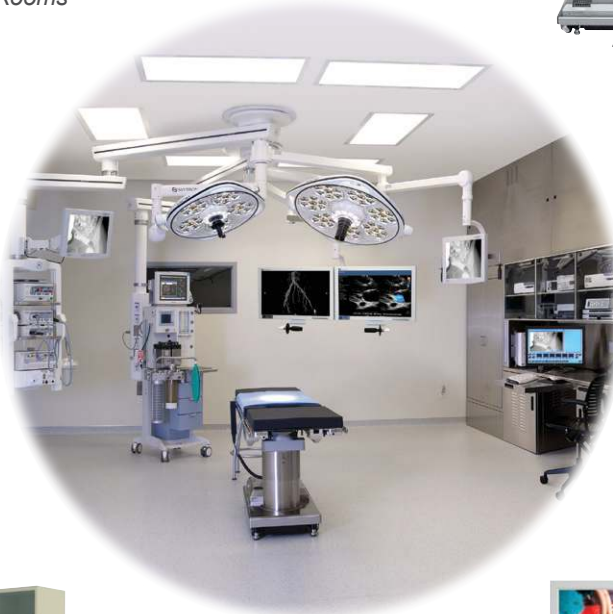

Space
Management
Booms



CSSD
Equipment



Fully Fiber Optic
Digital Integration



Solid Surface
Cabinets



Stainless



Asset Tracking and
Management Systems



HD Flat Panel Displays

This brochure contains information intended to illustrate a Skytron product. For questions about the availability of Skytron products in your area, contact your Skytron representative. All Skytron products are subject to change without prior notice.



5085 Corporate Exchange Blvd. SE • Grand Rapids, MI 49512
1-800-SKYTRON (759-8766) • email: sales@skytron.us

Local Representative

MAR-F-0002 REV4 12/13
PN L2-010-12